

Complete Mini / Mini Wide / Maxi Electrophoresis & Blotting System

Complete Mini / Mini Wide / Maxi Electrophoresis & Blotting System


The Complete Electrophoresis Systems include both modules for Gel Electrophoresis and Electroblotting and accessories, to provide a complete Mini, Mini Wide or Maxi gel casting, running and electroblotting system. The electroblotting module includes three inter-locking cassettes and six fiber pads.

A high intensity current is generated by coiled platinum electrodes. These in conjunction with the advanced cooling features allow rapid electroblotting in as little as one hour.


Features

- >Versatile: Interchangeable modules for slab gel and Electroblotting using a single universal buffer tank
- >Hinged Cassettes
- >High intensity electrodes
- >3 blot capacity
- >Rapid set-up cooling


MV-10CBS

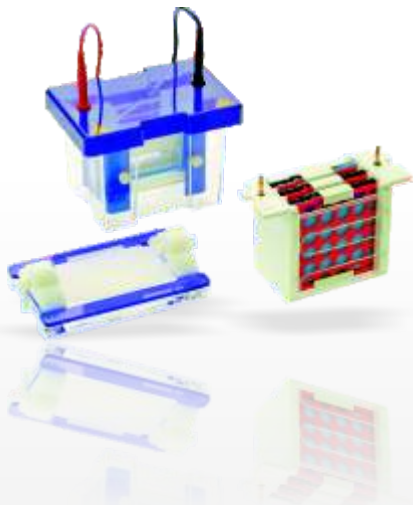

MV-10WCBS

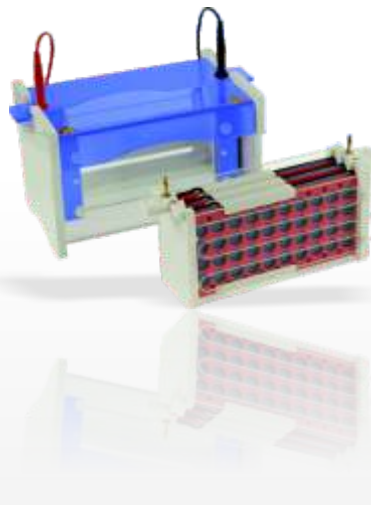


MV-20CBS

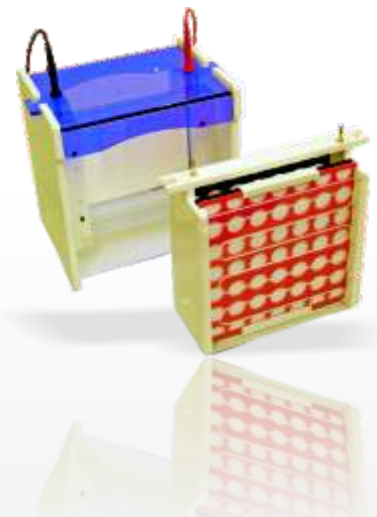

Ordering Information

MV-10CBS

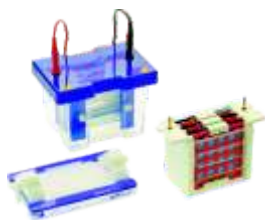


Complete system for Mini Vertical Electrophoresis & Blotting, comprising:

1 x Mini Vertical Unit, including
2 x 2mm thick notched glass plates
2 x 2mm thick plain plates with 1 mm thick bounded spacers, 1 x dummy plate
2 x combs (1 mm thick, 12 samples), 1 x casting base, silicone mat, cooling pack

1x Electroblotting module, including
1x Internal electroblotting module
3 x compression cassettes for gel sizes up to 100 x 100mm, 6 x fiber pads

MV-10WCBS

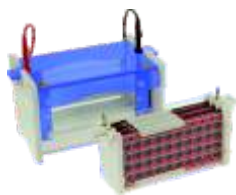

Complete system for Mini Wide Vertical Electrophoresis & Blotting, comprising:

1 x Mini Wide Vertical Unit, including
2 x 4mm thick notched glass plates
2 x 4mm thick plain plates with 1 mm thick bounded spacers, 1 x dummy plate
2 x combs (1 mm thick, 24 samples), 1 x casting base, silicone mat, cooling pack

1x Electroblotting module, including
1x Internal electroblotting module
3 x compression cassettes for gel sizes up to 200 x 100mm, 6 x fiber pads

MV-20CBS



Complete system for Maxi Vertical Electrophoresis & Blotting, comprising:

1 x Maxi Vertical Unit including
2 x 4mm thick notched glass plates
2 x 4mm thick plain plates with 1 mm thick bounded spacers, 1 x dummy plate
2 x combs (1 mm thick, 24 samples), 1 x casting base, silicone mat, cooling pack

1x Electroblotting module including
1x Internal electroblotting module
3 x compression cassettes for gel sizes up to 200 x 200mm, 6 x fiber pads



Contact Information

www.majorsci.com

info@majorsci.com

19959 Sea Gull Way
Saratoga, CA 95070
U.S.A

No.7, Hwai-Der Street,
Pan-Chiao City, Taipei Hsien 22045,
Taiwan

T/1-408-366-9866
F/1-408-446-1107

T/886-2-2256-6980
F/886-2-2256-6942

Innovative Life Sciences Tools